

Millipore®

Preparation, Separation,
Filtration & Monitoring Products

MILLIPORE SIGMA



FLASH® RT Allergen Indicator Protein Test


Cleaning Verification for Protein Removal and
Allergen Control Programs.

FLASH® RT is a total protein visual test that rapidly detects protein residues left on food contact surfaces after cleaning. Protein is a difficult food residue to remove. Most allergens are proteins, so quick verification of surface hygiene helps minimize the risk of cross-contact to allergen-free products. Using FLASH® RT regularly also helps reduce the opportunity for biological contamination.

As part of a HACCP allergen management program, FLASH® RT supports process verification requirements that ensure cleaning methods, independently validated to effectively remove allergens, are consistently applied.

Used together with the MVP ICON® for ATP testing, both can be incorporated into validated cleaning procedures (SSOPs) as part of a more comprehensive HACCP management program.

To find out more, visit:
SigmaAldrich.com

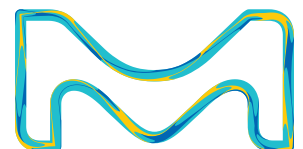

A Perfect Complement to ATP Testing

Colorimetric Assay. No Equipment Needed.

FLASH® RT is a self-contained sampling device containing a pre-moistened tip and all necessary reagents. Testing is done at ambient room temperature incubation and requires no equipment.

Manage Risks. Be Prepared. Everytime.


- Detects total protein including “Big 8” food allergens
- Rapid pass / fail result within 10 minutes
- Detects down to 20 µg total protein
- Simple online training available 24/7 with certification for audit compliance

3 Step Process. Simple.




SWAB

Sample a 4 x 4 inch
(10 x 10 cm) surface

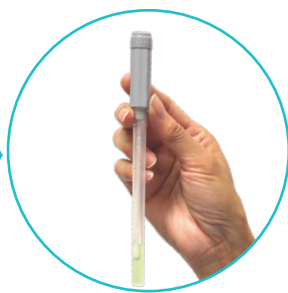

ACTIVATE

Activate swab and
shake for 5 seconds

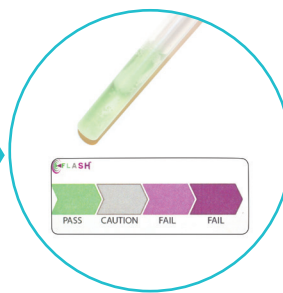

READ

Use chart to
interpret result
(Pass, Caution or Fail)

To place an order or receive technical assistance
In the U.S. and Canada, call toll-free 1-800-245-0113
For other countries across Europe and the world, please visit: [SigmaAldrich.com/offices](https://www.SigmaAldrich.com/offices)
For Technical Service, please visit: [SigmaAldrich.com/techservice](https://www.SigmaAldrich.com/techservice)

MilliporeSigma
400 Summit Drive
Burlington, MA 01803

[SigmaAldrich.com](https://www.SigmaAldrich.com)

© 2021 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. MilliporeSigma, the vibrant M, FLASH and MVP ICON are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

